

THE CYTOTOXICITY PATHWAY OF NATURAL KILLER CELLS IN CORD BLOOD COMPARED TO PERIPHERAL BLOOD

M. KHEIRANDISH,^{***} M. EBTEKAR,^{*} A.A. POURFATHOLLAH,^{***}
Z.M. HASSAN,^{*} AND A. KAZEM NEJAD^{*}

From the ^{*}Department of Immunology, School of Medical Sciences, Tarbiat Modarres University, Tehran,
and the ^{**}Iranian Blood Transfusion Research Center, Tehran, Iran.

ABSTRACT

The cytotoxic activity of natural killer cells is usually tested by radioactive assay (⁵¹Cr release assay), which detects the release of cytoplasmic contents after plasma membrane disintegration of dying cells. In contrast to this indirect evaluation of cytotoxicity, the assessment of cell damage by flow cytometry aims to provide a more exact characterization of the death pathway via detection of the percentage of apoptosis and necrotic cells. Annexin V-FITC (Axv-FITC) can be used to label cells in the early apoptotic state, while propidium iodide (PI) indicates late apoptosis or necrosis.

The NK cytotoxicity of cord blood (CB) and peripheral blood (PB) was determined after 4 hours of incubation in the absence of cytokines. After 4 hours *in vitro* incubation, co-staining with Annexin V-FITC (Axv-FITC) and propidium iodide (PI) permitted discrimination between viable, early apoptotic and necrotic cells. As we would expect, the cytotoxicity pathway in PB mononuclear cells (MNCs) consists of both apoptosis and necrosis pathways but in CB MNCs it almost consists of early apoptosis; and necrosis is negligible. With escalating E: T (effector: target) ratio changes in the percentage of apoptotic cells in PB samples were significantly higher than CB samples. The mechanism(s) of the low cytotoxicity of resting cord NK cells is not well understood. Complementary research in this field is recognized to elucidate the phenotypical and functional properties of CB cells and how they relate to maturational stages. CB studies are important for transplantation research and may provide insight to the suppressive mechanism by which the host-recipient could evade GVHD and rejection. *MJIRI, Vol. 18, No. 1, 55-60, 2004.*

Keywords: Peripheral Blood, Cord Blood, Natural Killer Cells, Apoptosis, Necrosis, Flow Cytometry.

INTRODUCTION

Human natural killer (NK) cells constitute a major population of lymphocytes distinguishable from B and T cells by a lack of antigen receptor. These cells play a

central role in the innate immune response in that they act both as accessory cells and as cytotoxic-effector cells, and as such, they are important in the surveillance and destruction of foreign, infected and tumor cells.¹

Unlike T and B cells, NK cells do not possess a specific antigen receptor with hypervariable region and kill target cells in a non-MHC restricted manner.²

NK cytotoxicity is mediated through two mechanisms:

Email: kheira_m2001@yahoo.co.uk

Cytotoxicity of NK Cells in Cord and Peripheral Blood

necrosis and apoptosis. Granules within the cytoplasm of NK cells contain at least two cytotoxic mediators: Perforin, which induces necrosis through membrane damage, and granzyme B, a serine esterase that splices nucleic acids.³

Direct induction of apoptosis by NK cells was controversial until recently, though it has now gained acceptance. Although molecules involved in signaling are not yet entirely defined, CD95 may be involved.⁴

NK cells also secrete soluble cytotoxic cytokines such as TNF α , which induces both necrosis and apoptosis *in vivo* and which may be more important in the function of NK cells than either of the direct mechanisms described above.⁵

Cord blood NK cells produce both perforin and granzyme B. The relative role of these mechanisms in cord NK - cell cytotoxicity is uncertain.^{6,7}

Cord blood NK cells, compared with adult peripheral blood NK cells, have a lower cytotoxic function against traditional cell lines such as K562 and Daudi.^{8,9}

However, the low cord blood NK cell cytotoxic function can be augmented by various cytokines, such as IFN- γ , IL-2, IL-7, IL-12 and IL-15 individually and in combination.^{8,10} IL-2, IL-12 and IL-15 raise cord blood NK cytotoxicity to the level of freshly isolated adult peripheral blood NK cells.¹¹

However, it remains lower than the level of adult cells treated with the same cytokines.^{8,9}

The cytotoxic activity of natural killer cells is usually tested by radioactive assay, which detects the release of cytoplasmic contents after plasma membrane disintegration of dying cells.

The ⁵¹Cr release assay is the most widely used method for measuring NK activity *in vitro* due to its reliability and simplicity.¹²

Nevertheless, as the usage of a radioactive material leads to various potential hazards, several groups have searched for easy and reliable alternatives.¹³

In contrast to this indirect evaluation of cytotoxicity, the assessment of cell damage by flow cytometry aims to provide a more exact characterization of the death pathway via detection of the percentage of apoptosis and dead cells.^{14,15}

One of the characteristic processes of the apoptosis pathway is a change in the plasma membrane architecture during early apoptosis; typical membrane compounds, such as phosphatidyl serine (PS) molecules, are redirected from the inner to the outer leaflet of the cell membrane without loss of membrane integrity.¹⁵ Annexin V-FITC (Axv-FITC), a molecule with high affinity for PS, can be used to label cells in the early apoptotic state, while propidium iodide (PI) indicates late apoptosis or necrosis.¹⁶

So far, we have described a fluorescence assay that

is based on the direct quantitative and qualitative flow cytometric analysis of cell damage at a single cell level. Assays were performed on 10 cord blood (CB) samples from full term normal deliveries. Results were compared with those obtained in PB samples from 10 normal adults.

MATERIAL AND METHODS

Effector cells

Isolation of mononuclear cells from peripheral and cord blood

Cord blood samples (N=10) were obtained from full term, normal delivery, healthy infants. Human peripheral blood samples (N=10) were obtained from healthy donors (18-45 years). Samples were collected in heparinized tubes and mixed with equal amounts of phosphate - buffer saline (PBS), pH=7.4.

Mononuclear cells (MNCS) were isolated by Ficoll-hypaque density - gradient centrifugation. Cell viability as measured by trypan blue exclusion was more than 99% for all cases.

Cells cultured in RPMI-1640 medium supplemented with 10% fetal calf serum (FCS, GIBCO, UK), 20 mmol/L, L-glutamine, 50 U/mL penicillin, 50 μ g/mL streptomycin and 20 mmol/L HEPES, referred to as complete medium. The cells were cultured at a density of 1×10^6 cells/mL in culture medium for the following assays.

Target cells

K562 cells (erythroleukemia cell line)(Iranian Blood Transfusion Organization, Cell Bank) as target cells were grown in complete medium as stationary cultures at 37°C in 5% CO₂. The cells were subcultured for 3 days before the assay, to be certain that they were in log phase.

NK Cytotoxicity assay

NK cell killing was assessed by modification of the method previously described.¹⁷ Target cells at the standard concentration and effector cells at various dilutions were added to create effector: target ratios (5:1, 10:1). Suspension cells were added to micro-plates and incubated for 4 hours; then, cultured cells were harvested and prepared for staining.

Staining with annexin-FITC and propidium iodide

After 4 hour co-incubation of target cells together with effector cells, total cell population was harvested from microtitre plates, collected into polystyrene tubes, then washed twice with PBS and high calcium - binding buffer (10 mM HEPES, NaOH, 150 mM NaCl, 1.8 mM CaCl₂, 1 mM MgCl₂, 5 mM KCl) then resuspended in 100 μ L binding buffer. Annexin-FITC (Axv-FITC, IQ-products, Netherlands) staining was performed with 5 μ L Axv-FITC and incubated for 10 min at room temperature in the dark.

Immediately before flow cytometric analysis 0.1µg PI (Sigma, USA) per 100µL of the binding buffer was added and single staining using Axv-FITC or PI alone was performed as controls.

FACS analysis

FACS analysis was carried out on a EPICS _XL Profile cytometer (Becton Dickinson, USA) by excitation at 488nm with an Argon laser, the emission of two fluorochromes was recorded through specific band pass filter: 530nm for FITC (FL1), > 675nm for PI (FL3). Specific fluorescence was analyzed after acquiring 5000 cells.

Target and effector gate and instrument variables were set according to forward and side scatter characteristics and using appropriate controls.

Analysis of fluorescence was performed after back

gating to exclude dead cells and aggregates.

Statistical analysis

Analyses of differences within and between groups were performed using the T-test for significance ($p < 0.05$).

RESULTS

Modeling for assaying apoptosis and necrosis

In order to distinguish target cells from effectors, the gating of target cells (K562) with forward and side scatter was performed. Data analysis was carried out by gating target cells followed by the analysis of Axv-FITC and PI-positive subpopulations. Axv-FITC single-positive cells (upper left quadrant) represented target cells in an early apoptotic state. Axv-FITC/PI double-positive

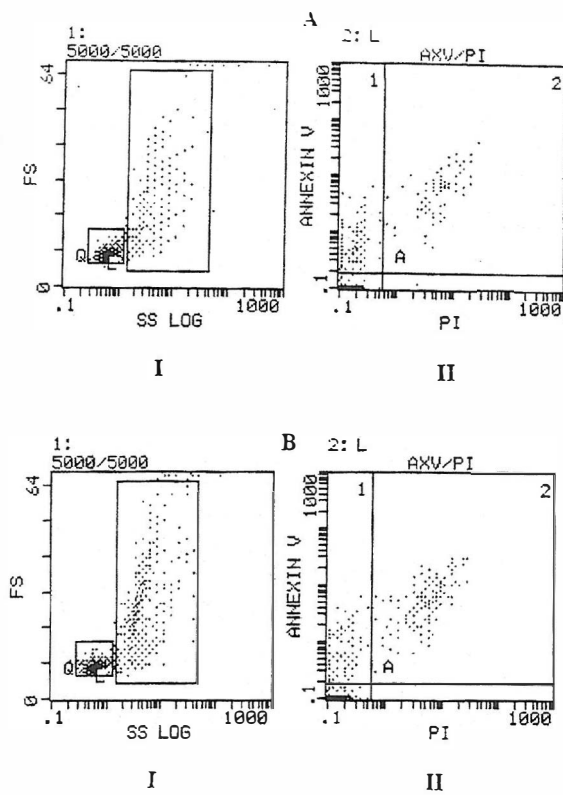


Fig. 1. Data acquisition. (I)-FCS/SS dot plot in an NK cytotoxicity assay with E: T ratios (A, 5:1 and B, 10:1) in peripheral blood. This dot plot representation discriminates between target (k562 cell line) and effector cells (peripheral blood mononuclear cells). (II)- A histogram plot of red fluorescence (x-axis: FL3, PI) and green fluorescence (Y-axis: FL1, annexin V-FITC). The percentage of FL-3 positive cells corresponds to the percentage of cell necrosis (lower right quadrant). The percentage of FL2-positive cells corresponds to the percentage of early apoptosis lysis (upper left quadrant). The percentages of FL2- and FL3- positive cells (dual positive) correspond to the percentage of late apoptosis and necrosis lysis (upper right quadrant).

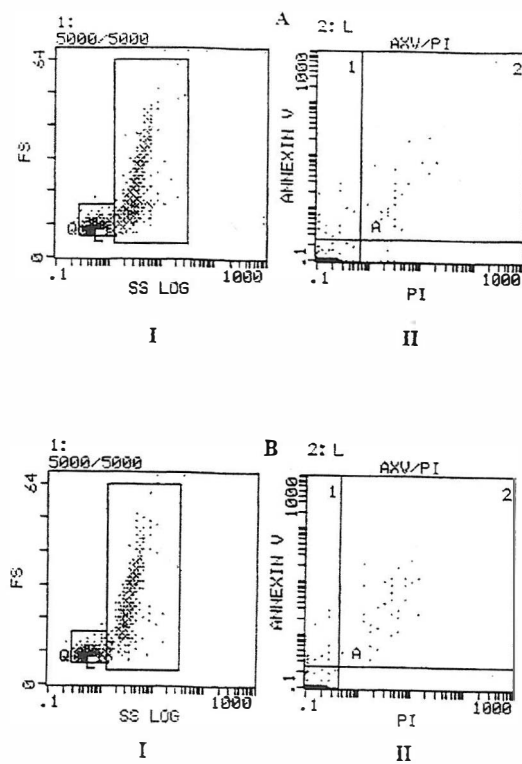


Fig. 2. Data acquisition. (I)-FCS/SS dot plot in an NK cytotoxicity assay with E: T ratios (A, 5:1 and B, 10:1) in cord blood. This dot plot representation discriminates between target (k562 cell line) and effector cells (cord blood mononuclear cells). (II)- A histogram plot of red fluorescence (x-axis: FL3, PI) and green fluorescence (Y-axis: FL1, annexin V-FITC). The percentage of FL-3 positive cells corresponds to the percentage of cell necrosis (lower right quadrant). The percentage of FL2-positive cells corresponds to the percentage of early apoptosis lysis (upper left quadrant). The percentages of FL2- and FL3- positive cells (dual positive) correspond to the percentage of late apoptosis and necrosis lysis (upper right quadrant).

Cytotoxicity of NK Cells in Cord and Peripheral Blood

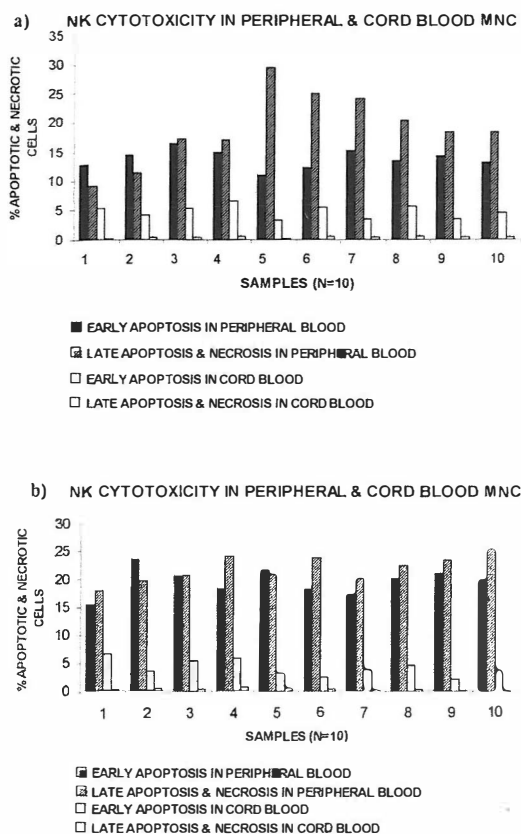


Fig. 3. Comparison NK cytotoxicity pathway of freshly isolated cord blood (CB) and peripheral blood (PB) mononuclear cells (MNCs). MNCs were isolated by Ficoll-hypaque density gradient centrifugation. Cells were washed and set up in a NK cytotoxicity assay at the effector: target (E: T) ratios (a) 5:1 and (b) 10:1. Data were compared using Student's t-test ($p < 0.05$).

tive cells (upper right quadrant) included cells in the late apoptotic and necrotic state. With escalating E: T ratios the percentage of Axv-FITC single positive or Axv-FITC double –positive cells is increased (Fig. 1, 2).

For the quantification of specific cytotoxicity, the population of non –specific positive cells (determined in the control sample without effector cells) in the corresponding quadrant has to be subtracted. The proportion of corrected Axv-FITC positive events (=Axv-FITC-positive event [sample] –Axv-FITC positive events [control]) to create total events (=total events [sample]–Axv-FITC positive events [control]) $\times 100$ reflects the real percentage of Axv-FITC positive cytotoxicity. The calculation was carried out separately for each quadrant and is exemplified for the Axv-FITC single – positive fraction.

As we would expect, the cytotoxicity pathway in PB

MNCS consisted of both apoptosis and necrosis pathways (in 5:1 ratios, the percentages are 13.9 ± 1.58 ; 19.14 ± 6.08 and in 10:1, 19.59 ± 2.33 ; 21.86 ± 5.63 ; respectively), as shown in Figure 3.

However, the cytotoxicity pathway in CB MNCS almost consists of early apoptosis and necrosis is negligible (in 5:1 ratios the percentages are 4.76 ± 2.22 ; 0.35 ± 0.16 ; in 10:1, 5.50 ± 1.64 and 0.32 ± 0.16 , respectively).

With escalating E: T ratios in CB NK, only the percentage of cell killing by early apoptosis pathway are increased.

Measurement of the early apoptosis induced by cord and peripheral natural killer cells

In order to estimate the early apoptosis induced by natural killer cells in both cord and peripheral blood, 10 samples were evaluated. The results in Table I show that a significant variation was noticed in cord blood compared to peripheral blood at ratios of 10:1 and 5:1. As shown in Fig. 3a, the percentage of cell killing by early apoptosis pathway in PB NK was significantly higher than in CB NK. As E: T ratios escalate, the percentage of apoptotic cells in PB samples significantly increases compared with CB samples.

Measurement of the late apoptosis and necrosis induced by cord and peripheral natural killer cells

In order to evaluate the late apoptosis and necrosis in natural killer cells in both cord and peripheral blood, 10 samples were evaluated, and the results show that a significant variation was noticed in peripheral blood comparing to cord blood at ratios of 10:1 and 5:1.

As shown in Fig. 3b, the percentage of cell killing by late apoptosis and necrosis pathway in PB NK was significantly higher than in CB NK. As E: T ratios escalate, the percentage of apoptotic cells in PB samples significantly increases compared with CB samples.

DISCUSSION

To date, over 1000 human umbilical CB transplants have been reported from different centers worldwide and it is generally agreed that CB represents an encouraging alternative to bone marrow (BM) transplantation.¹⁸

There are a variety of reasons for this, including the wider availability and easier access of CB as compared to BM. In addition, it has been suggested that there is a reduced graft –versus-host disease (GVHD) in CB transplants as compared to BM transplantation, even when donor and recipient are HLA-mismatched.^{19,22} The biological and immunological characteristics of cord blood cells differ from peripheral blood cells at several levels.

Cord blood lymphoid cells are characterized by naiveté, polyclonality of T cell repertoire, and absence

of certain molecules of activation, modification of cytokine production and receptors, and possible presence of suppressive cells. Whether these properties account for reduced capacity of transplanted cord blood cells to modulate a graft-versus-host disease remains to be determined.

It is clear that both T cells and NK cells have reduced function when isolated from CB as compared to adult cells and both of these cell types have been implicated in GVHD pathogenesis.^{23,26}

To support this, it has been documented that the NK cell phenotype within CB is different from those in the adult cells in that CB NK cells lack the CD57 antigen, which is found in 50-60% of adult NK cells.²³

However, whether this functional and phenotypic immaturity is sufficient to explain the apparently reduced GVHD seen in CB-transplanted individuals has not yet been answered.

Many studies indicate that cord blood NK activity is most likely suppressed by specific cells and others have also showed increased suppressor cells in cord blood.²⁷

The mechanism(s) of the low cytotoxicity of resting cord NK cells is not well understood.

There may be a mechanism intrinsic to cord blood NK cells, reflecting their immaturity, it may be due to extrinsic factors or it may be a combination of both.

Although CB NK cells subsets were deficient in or had less lytic activity against K562 cells compared to PB NK cells, they possess the lytic molecules perforin and granzyme B and when artificially stimulated to secrete their granules during lytic assays, were capable of lytic activity equivalent to that of PB NK cells.^{25,26}

A number of factors can modify the activity of NK cells to either increase or suppress cytotoxicity.²⁴

Umemoto et al. showed that in contrast to fresh resting cord blood NK cells, IL-2 activated NK cells induced both apoptosis and necrosis.²⁸

Complementary research in this field is recognized to elucidate the phenotypical and functional properties of CB cells and how they relate to maturational stages. CB studies are important for transplantation research and may provide insight in to suppressive mechanisms by which the host-recipient could evade GVHD and rejection.

Here, we have described a flow cytometry-based killing assay, which permits the fine evaluation of different stages of target cell death (early apoptosis and membrane damage) occurring during NK activity responses.

Target cells can be distinguished from effector cells using forward and sidescatter.

Co staining with Axv-FITC and PI permits the identification of target cells at different stages of cell death [early apoptosis (Axv-FITC+/PI-) and/ or necrosis (PI+)]. These results may provide further information on the death pathway of the target cells recognized by effector

cells. The assessment of cytolytic activity against target cells by Axv-FITC versus PI revealed different populations: Live cells (Axv-FITC-/PI), early apoptotic cells (Axv-FITC+/PI), late apoptotic/necrotic cells (Axv-FITC-/PI+ and/or Axv-FITC-/PI+).

Since Annexin V may be bound to phosphatidyl serine located either at the outer or the inner leaflet of the plasma membrane, the direct discrimination between late apoptotic and necrotic cells is critical.

With increasing E/T ratios, we observed an accumulation of two different populations of target cells (Fig. 1,2), the majority consisted of Axv-FITC+/PI- early apoptotic cells while the minority were found to be Axv-FITC-/PI+ late apoptotic target cells.

No significant increase in Axv-FITC-/PI+ necrotic cells was observed after 4 hours of incubation. These data suggest that most of the target cells recognized by effector cells die via the apoptosis pathway.

The flow cytometry assay has advantage over the ⁵¹Cr release assay; the assay may measure either apoptotic or necrotic cell death of the target cells.

Perforins are known to create channels in the cell membrane through which granzymes released by NK cells can induce the apoptotic cascade.²⁸ Therefore, some targets die a necrotic death from perforin lysis, while other cells die a slow, apoptotic death.

The ⁵¹Cr release assay measures the ⁵¹Cr released by necrotic cells, whereas apoptotic cells may not release ⁵¹Cr.

In conclusion, the flow cytometry assay is a quantitative method for the detection of cell-mediated cytotoxicity and does not use radioactive labeling. The reliability and reproducibility of this assay make it a useful alternative to the standard ⁵¹Cr release assay, providing a greater insight in to the mechanisms of cell-mediated cytotoxicity by different effector populations.

REFERENCES

1. Trinchieri G: Biology of natural killer cells. *Adv Immunol* 47: 187- 376, 1989.
2. Liao N, et al: Susceptibility to natural killer (NK) cells and impaired NK activity. *Science* 253: 199-202, 1991.
3. Smyth MJ, Trapani JA: Granzymes: exogenous proteinases that induce target cell apoptosis. *Immunol Today* 16: 202-6, 1995.
4. Vujanovic NL, Nagashima S, et al: Non-secretory killing by human NK cells. *J Immunol* 157: 1117-26, 1996.
5. Laster SM, Wood JG, et al: Tumor necrosis factor can induce both apoptotic and necrotic forms of cell lysis. *J Immunol* 141: 2629-34, 1997.
6. Berthou C, Legros-maidas S, et al: Cord blood T lymphocytes lack constitutive perforin expression in contrast to adult peripheral blood T lymphocytes. *Blood* 85: 1540-60, 1995.

Cytotoxicity of NK Cells in Cord and Peripheral Blood

7. Gaddy J, Broxmeyer HE: Cord blood CD16+D56- cells with low lytic activity are possible precursors of mature natural killer cells. *Cell Immunol* 180: 132-142,1997.
8. Lee SM, Suen Y, Chang L, et al: Decreased interleukin-12 (I1-12) from activated cord versus adult peripheral blood mononuclear cells and upregulation of interferon-gamma, natural killer, and lymphokine-activated killer activity by I1-12 in cord mononuclear cells. *Blood* 88: 945-54, 1996.
9. Qian JX, Lee Sm, Suen Y, et al: Decreased interleukin-15 from activated cord versus adult peripheral blood mononuclear cells and the effect of interleukin-15 in upregulating antitumor immune activity and cytokine production in cord blood. *Blood* 90: 316-17, 1997.
10. Chin TW, Ank BJ, et al: Cytotoxic studies in human newborns lessened allogeneic cell-induced (augmented) cytotoxicity but strong lymphokine-activated cytotoxicity of cord mononuclear cells. *Cell Immunol* 103: 241-51, 1996.
11. Ueno Y, Miyawaki T, Seki H, et al: Differential effects of recombinant human interferon-gamma and interleukin-2 on natural killer cell activity of peripheral blood in early human development. *J Immunol* 135: 180-4, 1985.
12. Brunner KT, Mauel J, Cerotini JC, et al: Quantitative assay of the cytolytic action of immune lymphoid cells on ⁵¹Cr-labeled allogeneic target cells: inhibition by isoantibody and by drugs. *Immunol* 14: 181-196, 1968.
13. Aubry JP, Blacke A, et al: Annexin used for measured apoptosis in the early events of cellular cytotoxicity. *Cytometry* 37: 197-204,1999.
14. Piriou L, Chilmonczyk S, et al: Design of a flow cytometric assay for the determination of natural killer and cytotoxic T-lymphocyte activity in human and in different animal species. *Cytometry* 41: 289-297, 2000.
15. Fischer K, Andreesen R, et al: An improved flow cytometry assay for the determination of cytotoxic T lymphocytes activity. *J Immunol Methods* 259: 159-164, 2002.
16. Martin SJ, Reutelings Perger CP, et al: Early redistribution of plasma membrane phosphatidyl serine is a general feature of apoptosis regardless of the initiating stimulus: inhibition by over expression of Bcl-2 and Abl. *J Exp Med* 1, 182(5): 1545-56, 1995.
17. Hatam L, Schuval S, et al: Flow cytometric analysis of natural killer cell function as a clinical assay. *Cytometry* 16: 59-66, 1994.
18. Hows JM: Status of umbilical cord blood transplantation in the year 2001. *J Clin Pathol* 54(6): 428-34,2001.
19. Vilmer E, Sterkers G, et al: HLA-mismatched cord blood transplantation in patient with advanced leukemia. *Bone Marrow Transplant* 2: 22,1991.
20. Flomenberg N, Keever CA: Cord blood transplants: potential utility and potential limitation. *Bone Marrow Transplant* 1: 115, 1992.
21. Miniero R, Ramenghi V: Umbilical cord blood stem cell transplantation. *Intern J Artificial Organs* 5: 113, 1993.
22. Keever C, Abu-Hajir M, et al: Characterization of the alloreactivity and anti-leukemia reactivity of cord blood mononuclear cells. *Bone Marrow Transplant* 15: 407, 1995.
23. Cohen SBA, Dominiguez E, et al: The immunological properties of cord blood: overview of current research Presented at the 2nd Eurocord Workshop. *Bone Marrow Transplantation* 22, (1): s22-s25, 1998.
24. Domingues E, Madrigal JA: Fetal natural killer cell function is suppressed, *Immunology* 94: 109-114, 1998.
25. Gardiner CM, O'meara A: Differential cytotoxicity of cord blood and bone marrow derived natural killer cells. *Blood* 91: 1, 207-213, 1998.
26. Han P, Dominguez E, et al: Immunobiology of neonatal natural killer cells in: Cohen SBA, Gluckman E, et al. (eds.), *Cord blood Characteristics: Role in stem cell transplantation*. London: Martin Dunitz, pp. 61-72, 2000.
27. El marsafy S, Dosquet C, et al: Study of cord blood natural killer suppressor activity. *Eur J Haematol* 66(4): 215-20, 2001.
28. Umemoto M, Azuma E, et al: Two cytotoxic pathways of natural killer cells in human cord blood implications in cord blood transplantation. *Br J Haematol* 98: 1037-40, 1997.